

HPV Self Sample Collection Kit

Take Control of Your Cervical Health - Privately, Easily, at Home.

Why HPV testing?

99.7% of cervical cancer cases are linked to persistent high-risk HPV infection



Cervical cancer is preventable – through regular screening for early detection and treatment



WHO key recommendations^[1] on HPV testing and self-sampling^[2]:

WHO Guidance	Key Message
✓ HPV DNA test as primary screening	More sensitive, allows 5–10 year screening intervals
✓ Self-sampling is reliable	WHO endorses vaginal self-sampling for effective cervical cancer screening
✓ Start screening at age 30	Every 5–10 years with HPV test
✓ Self-sampling increases coverage	Especially for women with limited access to screening services

When to use this kit



Living far from screening centers



Uncomfortable with pelvic exams or concerned about privacy



Need long-term or long-distance HPV sample preservation

Sources:

[1]World Health Organization. (2024). WHO guideline for screening and treatment of cervical pre-cancer lesions for cervical cancer prevention: use of dual-stain cytology to triage women after a positive test for human papillomavirus (HPV). 2nd ed. Geneva: World Health Organization.

[2]World Health Organization. (2024). WHO Information Brochure on Self-Care Interventions for Human Papillomavirus (HPV) self-sampling as part of cervical cancer screening and treatment.




Kit Components

- Transport Medium - 1
- Gynecological sampler - 1
- ID Label - 1
- Biosafety Bag - 1
- Instruction for Use - 1



Sampling Operations

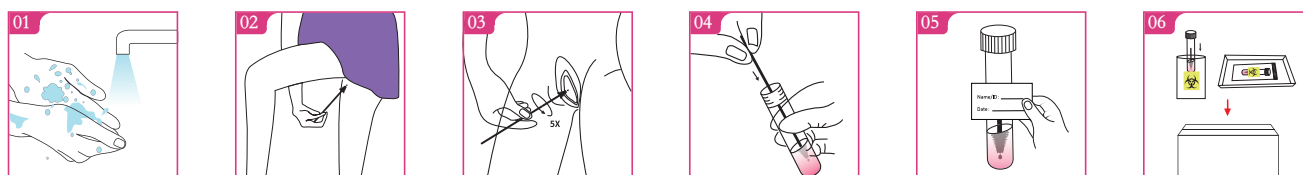

★ Single use only. Do not reuse or share. Results provided by your local medical laboratory.

Product Features

- 01 The preservation solution contains antibiotics to stabilize virus and prevent bacteria/fungi.
- 02 Works with standard HPV PCR lab tests.
- 03 Storage & transport summary

Condition	Maximum duration
Room temperature (2-30°C)	7 days (shipping)
Refrigerated (2-8°C)	15 days
Frozen (-20±5°C)	6 months

- Avoid: extreme heat, repeated freeze-thaw cycles

Ordering Information

Cat. No.	Product Name	Kit Size	Certificate
PH03010101T2001	Human Papillomavirus Self Sample Collection Kit	1 Set/Kit	CE-IVDR

- Expiration | 2 years
- Storage Condition | 2-30°C

Ver: 0.1 Date: 2026-05-25



Jiangsu Mole Bioscience Co., Ltd.
No.1, Lingzhi Road, Medical New & Hi-tech
District (Gaogang District), 225300 Taizhou,
Jiangsu Province, China

Tel: 0571-87209310
Email: info@molechina.com
www.molechina.com